

LITEAID

Neck Massager LA-045P



Instruction Manual

New York, NY

This manual should be kept on hand for easy reference. It is your guide to a safe and efficient operation.

Customerservice@liteaid.com
www.liteaid.com

Product: Orion Elite Neck Massager - Pad
Model Number: LA-045P

Weight: 630 g
Power: Instrument 6VDC≤6W
Controller CR2032 Button Battery 3VDC
Electric Impulse Frequency: 1- 1000 Hz
Width of Electric Impulse: 20-400

MADE IN CHINA

I. Warnings

This product is intended only for use by adults. It is not a toy.

1. Please see a physician before use if you are pregnant or if you suffer from severe neck or back pain, or from other prior medical conditions.
2. Always turn off your neck massager when not in use.
3. Do not use your neck massager near water sources such as sinks, bathrooms, or showers.
4. Store neck massager in a safe and dry environment.
5. Do not expose your neck massager to extreme cold or extreme heat.
6. Never use units that have been broken or damaged.
7. Always practice caution when using massaging devices. If you feel dizzy or out of breath, always consult with your physician prior to use.
8. Cannot be used with a PACEMAKER
9. Do not use near a microwave or other infrared emitting devices.
10. Not recommended for use by individuals with skin sensitivity or heat sensitivity.
11. Do not use when sleeping or showering.
12. Do not attempt to repair if broken.

Wipe clean with a soft cloth.

Store in a dry temperature controlled area.

II. Introduction:

Fig. 2


This neck massage therapy was adapted from Chinese Meridian beliefs. Designed to conform to the curve of the neck, this device combines low frequency electric therapy, magnetic therapy, far infrared therapy and vibration massage therapy.

III. Features:

- Designed specifically to conform to the curve of the neck without discomfort.
- Far infrared heating, electric impulse, and vibration massage.
- Temperature settings according to personal need and preference.
- 6 Vibration modes
- Wireless remote and backlit LCD screen.

IV. Product Diagram

Fig. 1

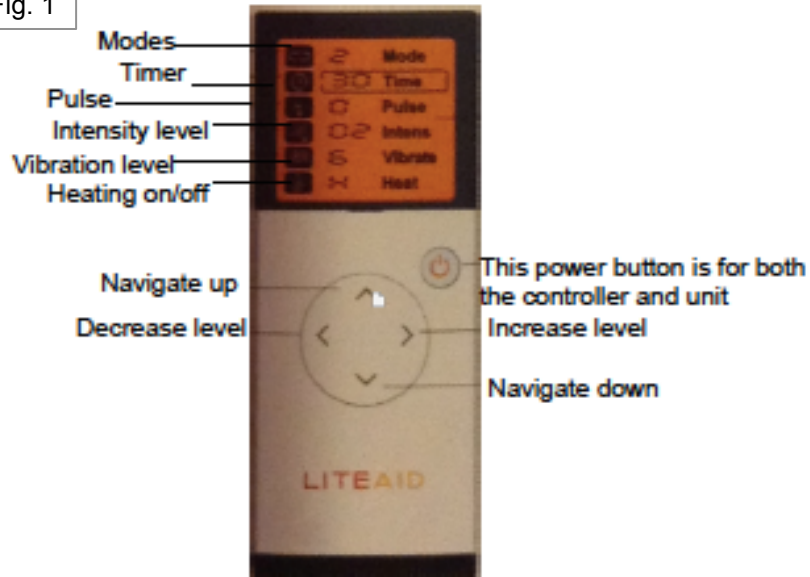

Steady light: ready for use

Rapid flashing light: charging


Slow flash: fully charged

LED Power Light



Fig.3


Infrared receiver

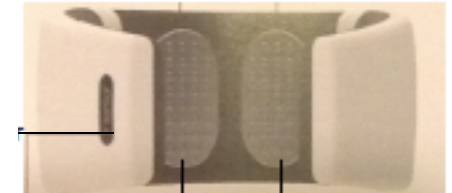

Heating and Electronic Impulse Plate

Steady light:
ready for use

Fig.4

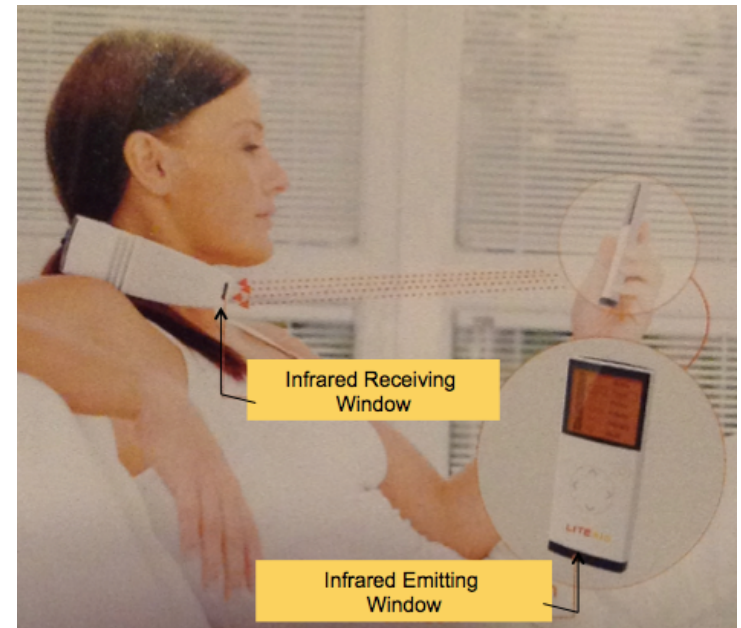

Important: In order to use the neck massager, the remote must be facing (see below) the transmitter and within 2 feet of the neck massager while changing functions using the remote.

V. Before Setup

1. At initial setup, keep the infrared emission window on the controller and the infrared light receiving window on the main unit directly opposite each other, as shown in Figure 4. This position is necessary for the remote controller to operate. Each valid operation of the remote will beep when starting.
2. Charging: The unit must be charged with regular usage. To charge, switch the ON/CHARGE button to "CHARGE," and the indicator will flash green. When the indicator light stops flashing and remains a steady green light, the unit is fully charged. It will take between 3-4 hours to fully charge the neck massager.

VI. Attaching the EMS Pads:


Each neck massager comes with 1 pairs of connectivity silicon pads. These pads are designed to gently conduct the electronic pulse for a comfortable user experience.

- Each set of pads is suitable for 3-4 months of usage.

To attach:


1. Remove silicon pads from plastic bag.
2. Carefully remove plastic liner on the side A of pads (see image below)

Fig. 5

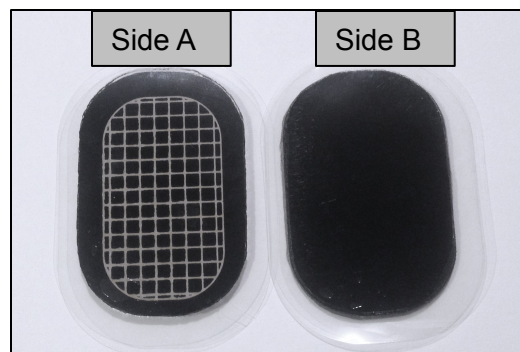


3. Being careful not to touch adhesive surface, place the pad on the heating plate, leaving side B facing outward with plastic liner attached.
4. Repeat steps 2 and 3, and place plastic liners back in plastic bag for further use.
5. Before use, remove plastic liners on side B and place neck massager on clean and dry neck. Never use the massage pads or neck massager if you have open sores or cuts on the neck.

VII. To Operate:

1. Set the operation switch to the "ON" position
2. Place neck massager on neck (read Section VI prior to doing so).
3. Press the  button to begin operation. Note the buzz and light indicating start-up. The LCD on the remote will indicate the mode while the massager is in use. When the remote and the massager synchronize, the unit will buzz.
4. Use the \wedge and \vee buttons to switch between modes,
5. Use the $>$ button to increase and the $<$ button to decrease.

VIII. MODES:

1. To choose modes, first use the \wedge and \vee in order to select 
2. To select between the 6 modes, press the $>$ and $<$ buttons.

Mode 1: AUTO MODE

- You can adjust the intensity of vibration, and the type of vibration pattern (intensity, and vibrate respectively). If you let this mode continue, it will take you through all vibration patterns and pulse patterns, allowing you to adjust the intensity as you need.

Mode 2: HEAT AND VIBRATION

- Maximum intensity is 3 with this vibration mode, with heat on high intensity. Adjustable vibration only – heat will remain the same throughout.

Mode 3: PULSE AND HEAT

- Adjust both heating intensity and pulse intensity. No vibration.

Mode 4 : HEAT

- Mode is heat only. Adjust H (H is high, L is low, 0 is off)

Mode 5: VIBRATION

- Vibration type and intensity can be adjusted. Heat is OFF.



Mode 6: PULSE

- Pulse with NO heat. You can adjust intensity as necessary, but no vibration or heat.

IX. Basic Commands

1. Intensity – Vibration and Electrical Impulse Manual Settings



1. Use the \wedge and \vee to choose the intensity symbol.
2. Press the $>$ and $<$ buttons in order to increase the intensity on the vibration (3 levels of intensity) or the electrical impulse
3. Both Electrical impulse and Vibration can be used manually with with Intensity button, or on specified modes with their individual mode settings. See section 3 for vibration and section 4  for electric pulse modes. 

★ IMPORTANT:

Always use caution when increasing the electrical impulse feature. The intensity can be increased up to level 16. It is recommended that new users to start at 1 and test in increments of one until a comfortable level is reached. Such levels are generally between 5 and 6.

★ Please note that to rapidly lower intervals the remote needs to be pointing at the infrared emission window (see Figure 4) and will beep with each change and will take a few seconds for each change to be accomplished.

2. Increasing and Decreasing Heat



1. Use the \wedge and \vee to choose the heat symbol.
2. Press the $>$ and $<$ in order to increase the heat to high (H), low (L) or none (0).

3. Selecting Vibration Modes



1. Use the \wedge and \vee to choose the vibration symbol.
2. Press the $>$ and $<$ in order to select between the 6 vibration modes:
 1. Auto-vibration
 2. Vibration only
 3. Knead
 4. Beating
 5. Press
 6. Rolling

4. Electric Impulse Modes



1. Use the \wedge and \vee to choose the electrical impulse symbol.
2. Press the $>$ and $<$ in order to select between the 6 impulse modes.

5. Timer



1. Use the \wedge and \vee to choose the timer symbol.
2. Press the $>$ and $<$ in order to set the timer within 5 minute intervals.

★ It is not recommended using the neck massager for longer than 15 minutes at a time.

X. Troubleshooting

Problem	Potential Solution
Neck massager does not turn on	Leave neck massager on charger for at least 3-4 hours
Neck massager will not switch modes or switches to rapidly.	Allow 15 seconds after increasing/decreasing intensity or changing mode in order to feel the change in intensity.
Neck massager will not connect with remote	See Figure 4. Check that infrared emitting window on remote is close to infrared receiving window on neck piece.

★ IMPORTANT:

Always use caution when increasing the electrical impulse feature. The intensity can be increased up to level 16. It is recommended that new users to start at 1 and test in increments of one until a comfortable level is reached. Such levels are generally between 5 and 6.

★ Please note that to rapidly lower intervals the remote needs to be pointing at the infrared emission window (see Figure 4) and will beep with each change and will take a few seconds for each change to be accomplished.

Order Replacement Connectivity silicon pads

<http://liteaid.com/products/replacement-pads-of-orion-elite-neck-massager>

Questions about your Liteaid USA Product?

Contact us at:
(646) 504 43 94
customerservice@liteaid.com